

Product Catalog.

Dental resins, equipment & accessories.



TCS

An ISO 13485 Certified Company

TCS is a dental manufacturing company dedicated to providing dental laboratories with exceptional quality materials, equipment, and hands-on training courses that are accessible and reliable. At the core of our brand, we believe that one of the key components to creating successful removable appliances is through forming strong partnerships, because for TCS, manufacturing is more than just a business, it is a collectivist culture.

2000

FOUNDED IN CA, USA

3

GENERATIONS

70+

COUNTRIES WE SELL TO



The TCS Experience

Our partnership at a glance

Dynamism. An energetic and passionate team moves forward together and constantly progresses to add value to our customers.

Education. Outstanding technical support and hands-on training designed to grow confidence in creating exceptional TCS dental appliances.

Innovation. A philosophy to continuous improvement guides us to develop better methods, tools, materials, educational programs, etc.

Manufacturing In-House. The majority of our products, including our resins and equipment are manufactured in-house in CA, USA.

Together We Smile - TCS Charitable Events. Our industry is fortunate to be part of a community where giving back literally means putting smiles on people's faces.

Win - Win Philosophy. Your success is our success. We aim for long term relationships based on mutual trust, respect and benefit.



Table of Contents

Unbreakable™	6
iFlex™	10
Karadent™	14
Digital Discs.....	18
Injection Units.....	20
Basic Kits.....	21
Complete Kits.....	22
JP90.....	26
Digital Furnace.....	28
Rental Equipment.....	29
Waxes / Flasks.....	30

Spatulas / Knives.....	31
Processing Tools.....	32
Finishing Tools.....	33
Polishing Tools.....	34
TCS Fresh.....	36
Get More With TCS.....	38
Technical Support.....	40
Online Training.....	41
Promote Your Department.....	42
Together We Smile.....	44

Unbreakable™


Flexible Partial



Monomer Free
Hypoallergenic



- Extremely resilient, high performance nylon thermoplastic.
- Guaranteed unbreakable due to its unique properties.
- Designed for exceptional memory to hold its form.
- Offers an ideal degree of flexibility.
- A modern fit, lightweight and durable.

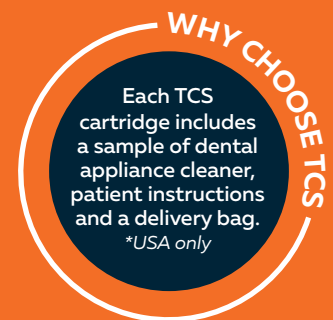

Unbreakable™

TCS resins are available in bulk or aluminum cartridges that come in various sizes/diameters.

Be sure to order the correct diameter cartridge to fit your injection unit.

Refer to [youtube.com/tcsdentalinc](https://www.youtube.com/tcsdentalinc) (TCS Short Videos Playlist).

Millable discs available for CAD/CAM - Refer to pages 18-19.



Shade	Size	25 mm	28 mm	22 mm
Light Pink (Shade 1)	S	4002-01/S1		
	M	4000-01/M1	4000-01/MF1	4005-01/M22-1
	L	4005-01/L1	4005-01/LF1	4005-01/L22-1
Standard Pink (Shade 2)	S	4002-01/S2		
	M	4000-01/M2	4000-01/MF2	4005-01/M22-2
	L	4005-01/L2	4005-01/LF2	4005-01/L22-2
Light / Dark Pink (Shade 3)	S	4002-01/S3		
	M	4000-01/M3	4000-01/MF3	4005-01/M22-3
	L	4005-01/L3	4005-01/LF3	4005-01/L22-3
Dark Pink (Shade 4)	S	4002-01/S4		
	M	4000-01/M4	4000-01/MF4	4005-01/M22-4
	L	4005-01/L4	4005-01/LF4	4005-01/L22-4
Natural (Shade N)	S	4002-01/SN		
	M	4000-01/MN	4000-01/MFN	4005-01/M22-N
	L	4005-01/LN	4005-01/LFN	4005-01/L22-N


Each pack includes the following:

5 delivery bags | 5 patient care instructions | 5 samples of TCS Dental Appliance Cleaner (1-week supply).*

Each kilo includes the following:

50 delivery bags, 50 patient care instructions.

*Not applicable to JP Bullets. USA only.**

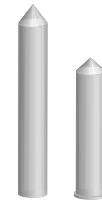

25 mm w / tab	25 mm w / nozzle	1 Kilo	500 Grm	JP Bullets
				4002-JP/S1
4005-01/MT25-1		4003-01B/1	4003-01B/1-.5	
4005-01/LT25-1	4005-01/LNZ25-1			4005-JP/L1
				4002-JP/S2
4005-01/MT25-2		4003-01B/2	4003-01B/2-.5	
4005-01/LT25-2	4005-01/LNZ25-2			4005-JP/L2
				4002-JP/S3
4005-01/MT25-3		4003-01B/3	4003-01B/3-5	
4005-01/LT25-3	4005-01/LNZ25-3			4005-JP/L3
				4002-JP/S4
4005-01/MT25-4		4003-01B/4	4003-01B/4-.5	
4005-01/LT25-4	4005-01/LNZ25-4			4005-JP/L4
				4002-JP/SN
4005-01/MT25-N		4003-01/BN	4003-01/BN-.5	
4005-01/LT25-N	4005-01/LNZ25-N			4005-JP/LN


iFlex[™]
Flexible Partial



Monomer Free
Hypoallergenic



- Highly translucent, polyolefin thermoplastic.
- Simple to finish and polish.
- Stain resistant and easy to maintain.
- Clinically unbreakable.
- Offers an ideal degree of flexibility.
- A modern fit, lightweight and durable.

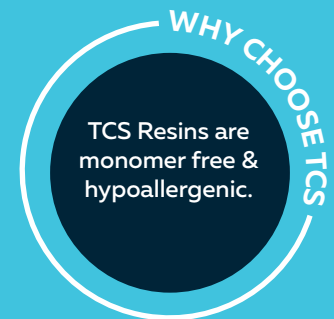

iFlex™

TCS resins are available in bulk or aluminum cartridges that come in various sizes/diameters.

Be sure to order the correct diameter cartridge to fit your injection unit.

Refer to [youtube.com/tcsdentalinc](https://www.youtube.com/tcsdentalinc) (TCS Short Videos Playlist).

Millable discs available for CAD/CAM - Refer to pages 18-19.



Shade	Size	25 mm	28 mm	22 mm
Light Pink (Shade 1)	S	5002-01/S1		
	M	5000-01/M1	5000-01/MF1	5005-01/M22-1
	L	5005-01/L1	5005-01/LF1	5005-01/L22-1
Standard Pink (Shade 2)	S	5002-01/S2		
	M	5000-01/M2	5000-01/MF2	5005-01/M22-2
	L	5005-01/L2	5005-01/LF2	5005-01/L22-2
Light / Dark Pink (Shade 3)	S	5002-01/S3		
	M	5000-01/M3	5000-01/MF3	5005-01/M22-3
	L	5005-01/L3	5005-01/LF3	5005-01/L22-3
Dark Pink (Shade 4)	S	5002-01/S4		
	M	5000-01/M4	5000-01/MF4	5005-01/M22-4
	L	5005-01/L4	5005-01/LF4	5005-01/L22-4
Natural (Shade N)	S	5002-01/SN		
	M	5000-01/MN	5000-01/MFN	5005-01/M22-N
	L	5005-01/LN	5005-01/LFN	5005-01/L22-N



Each pack includes the following:

5 delivery bags | 5 patient care instructions | 5 samples of TCS Dental Appliance Cleaner (1-week supply).*

Each kilo includes the following:

50 delivery bags, 50 patient care instructions.

*Not applicable to JP Bullets. USA only.**

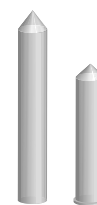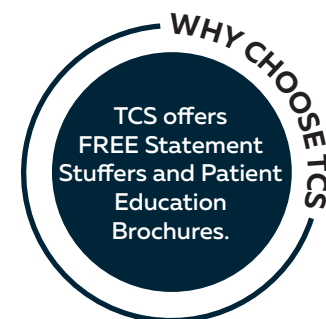

25 mm w / tab	25 mm w / nozzle	1 Kilo	500 Grm	JP Bullets
--				5005-JP/S1
5005-01/MT25-1		5003-01B/1	5003-01B/1-.5	--
5005-01/LT25-1	5005-01/LNZ25-1			5005-JP/L1
--				5005-JP/S2
5005-01/MT25-2		5003-01B/2	5003-01B/2-.5	--
5005-01/LT25-2	5005-01/LNZ25-2			5005-JP/L2
--				5005-JP/S3
5005-01/MT25-3		5003-01B/3	5003-01B/3-.5	--
5005-01/LT25-3	5005-01/LNZ25-3			5005-JP/L3
--				5005-JP/S4
5005-01/MT25-4		5003-01B/4	5003-01B/4-.5	--
5005-01/LT25-4	5005-01/LNZ25-4			5005-JP/L4
--				5005-JP/SN
5005-01/MT25-N		5003-01B/N	5003-01/BN-.5	--
5005-01/LT25-N	5005-01/LNZ25-N			5005-JP/LN

Karadent™


Full Dentures
Clear Frameworks



Monomer Free
Hypoallergenic



- Semi-rigid, microcrystalline polymer thermoplastic.
- Ideal for full dentures, semi-flexible partial dentures, clear clasps, frameworks, and nightguards.
- An upgrade from traditional acrylic dentures yet as simple to finish and polish.
- Finishes completely transparent in our Clear shade.
- A modern fit, lightweight and durable.

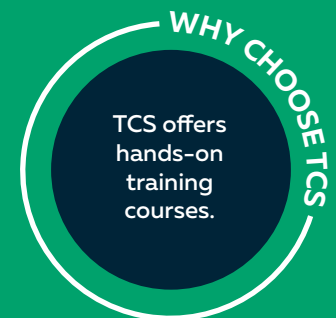


Karadent™

TCS resins are available in bulk or aluminum cartridges that come in various sizes/diameters.

Be sure to order the correct diameter cartridge to fit your injection unit.

Refer to [youtube.com/tcsdentalinc](https://www.youtube.com/tcsdentalinc) (TCS Short Videos Playlist).

Millable discs available for CAD/CAM - Refer to pages 18-19.

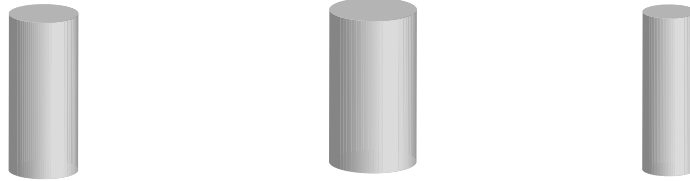

Shade	Size	25 mm	28 mm	22 mm
Standard Pink (Shade 2)	S	2000-01/S2	--	--
	M	2000-01/M2	2000-01/MF2	2005-01/M22-2
	L	2005-01/L2	2005-01/LF2	2005-01/L22-2
Dark Pink (Shade 4)	S	2000-01/S4	--	--
	M	2000-01/M4	2000-01/MF4	2005-01/M22-4
	L	2005-01/L4	2005-01/LF4	2005-01/L22-4
Natural (clear) (Shade N)	S	2000-01/SN	--	--
	M	2000-01/MN	2000-01/MFN	2005-01/M22-N
	L	2005-01/LN	2005-01/LFN	2005-01/L22-N



Each pack includes the following:

5 delivery bags | 5 patient care instructions | 5 samples of
TCS Dental Appliance Cleaner (1-week supply).*

Each kilo includes the following:

50 delivery bags, 50 patient care instructions

*Not applicable to JP Bullets. USA only.**

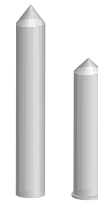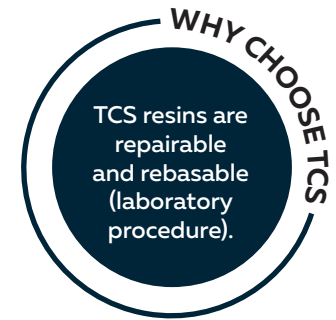

25 mm w / tab	25 mm w / nozzle	1 Kilo	500 Grm.	JP Bullets
--	--			2002-JP/S2
--	--	2003-01/B2	2003-01B/2-.5	--
--	2005-01/LNZ25-2			--
--	--			2002-JP/S4
--	--	2003-01/B4	2003-01B/4-.5	--
--	2005-01/LNZ25-4			--
--	--			2002-JP/SN
--	--	2003-01/BN	2003-01/BN-.5	--
--	2005-01/LNZ25-N			--


Dive Into the Digital Age



Available Soon



Millable Discs *for a modern fit*

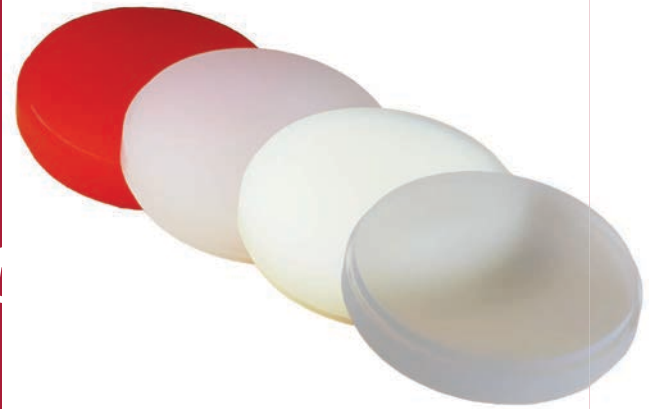

- Digital technology allows for ultimate precision.
- Appliances come out with a smooth texture, close to polishing.
- Polishing process is quick (approx. 2 minutes).
- Milling an appliance from a disc is kinder on your tools.

Contact us for details
(866) 426-2970

TCS Injection Units.

TCS offers a variety of Injection Units to fit laboratories preferences and budget.

MOST POPULAR



Hydraulic Injector (electric).

Equal pressure with every injection.

Minimal maintenance.

No air connections required.

Requires an external heater
(TCS Digital Furnace)

Injector only 110V 3504-01

Injector only 220V 3504-02

SPECS.

Yields approx. 30 cases per day
(up to 120 cases with additional
furnaces).

H: 22" **L:** 11" **W:** 6.5"

50 lbs. (23 kilos).

2 year warranty
Made in USA

MOST ECONOMICAL



Manual Injector.

Minimal maintenance.

Built to last.

Simple to use.

Requires an external heater
(TCS Digital Furnace)

Injector only 110V 3506-01

SPECS.

Yields approx. 30 cases per day
(up to 120 cases with additional
furnaces).

H: 22" **L:** 22" **W:** 17"

42 lbs. (19 kilos).

2 year warranty
Made in USA

ALL IN ONE



Automatic Injector.

Compact design with a built-in furnace.

Equal pressure with every injection.

Injects with a touch of a button.

Requires hookup to air pressure (min 140 psi)

Injector only 110V 3525-01

Injector only 220V 3525-02

SPECS.

Yields approx. 30 cases per day.

H: 23" **L:** 8" **W:** 13.5"

40 lbs. (18 kilos).

1 year warranty
Made in USA


Basic Kits.

Kits include the basic essentials to inject TCS Resins.

Hydraulic Basic Kits.

110V 3511-01

220V 3511-02

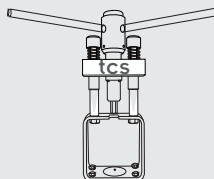

Manual Basic Kits.

110V 3508-01

220V 3508-02



1 digital furnace



1 cylinder with handle



1 knock out rod
1 knock out stand



1 thick bronze plug
1 thin bronze plug




1 injection flask



2 cartridges of each:
Unbreakable™, iFlex™,
Karadent™



1 mold release spray

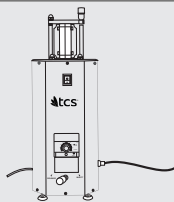

Automatic Air Basic Kits.

110V 3528-01

220V 3528-02



1 injection flask



2 cartridges of each:
Unbreakable™, iFlex™



1 mold release spray

Note: It is not recommended to inject Karadent™ by tcs® using the Automatic Air Injector.



TCS Complete Injection Kits.

Includes all the essentials necessary to process TCS Resins.

Hydraulic Complete Kits.

110V 3501-01

220V 3501-02

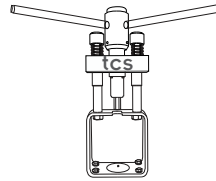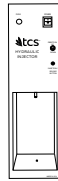

Manual Complete Kits.

110V 3505-01

220V 3505-02



2 pks - Unbreakable™
2 pks - iFlex™
2 pks - Karadent™



1 digital furnace



2 cylinders with handle



1 knock out rod
1 knock out stand



2 thick bronze plug
1 thin bronze plug



2 injection flasks



10 pcs. palate wax



10 pcs. lingual wax



10 pcs. sprue wax



3 pcs. 9.5 sprue wax



1 carving knife



1 cement



1 EZ flow sample



1 mold release spray



1 fusing liquid



1 hot air welder



Twist drill samples



TCS putty (400 gms)

1 econo cutter
1 econo mandrel1 white knock down
wheel & chuck1 rubberized abrasive 304
1 rubberized abrasive 306

10 pink grinding wheels



10 brown rubber wheels



5 green mounted stones

2 large rubber points
2 small rubber points

2 mini fuzz buffs



1 pink zirconia coarse bur



1 lb polycril



1 turbo buff



4 rag buffs


sample pieces of
polishing paste1 sample model
upon request

Items in the kits may change without notice. We strive for continuous improvement, if we discover a product that works better with our resins, we may switch out an item(s).



TCS Complete Injection Kits.

Includes all the essentials necessary to process TCS Resins.

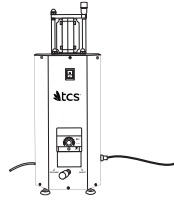

Automatic Complete Kits.

110V 3526-01

220V 3526-02



3 pks - Unbreakable™
3 pks - iFlex™



2 injection flasks



10 pcs. palate wax



10 pcs. lingual wax



10 pcs. sprue wax



3 pcs. 9.5 sprue wax



1 carving knife



1 cement



1 EZ flow sample



1 mold release spray



1 fusing liquid



1 hot air welder



Twist drill samples



TCS putty (400 gms)



1 econo cutter
1 econo mandrel



1 white knock down
wheel & chuck



1 rubberized abrasive 304
1 rubberized abrasive 306



10 pink grinding wheels



10 brown rubber wheels



5 green mounted stones



2 large rubber points
2 small rubber points



2 mini fuzzi buffs



1 pink zirconia coarse bur



1 lb polycril



1 turbo buff



4 rag buffs




sample pieces of
polishing paste



1 sample model
upon request

Note: It is not recommended to inject Karadent™ by tcs® using the Automatic Air Injector.

Items in the kits may change without notice. We strive for continuous improvement, if we discover a product that works better with our resins, we may switch out an item.

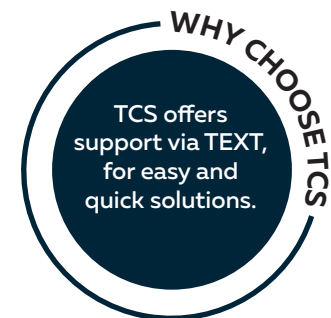

JP90 Handheld Injector.



Used in conjunction with the TCS Digital Furnace. Sold separately.

Injects:

- » 1-3 tooth unilaterals.
- » Up to 5 consecutive tooth partials/saddles, saddles for combination cases.
- » Repairs.

JP90 Injector only

7000-JP

SPECS.

H: 6" L: 2-1/4" W: 11"

2 lbs. (1 kilo).

1 year warranty

Made in USA



YouTube: [youtube.com/tcsdentalinc](https://www.youtube.com/tcsdentalinc)
(TCS Equipment playlist)



JP Flask.

Fabricated from high precision stainless steel to withstand warpage under extreme heat.

7005-JP



Metal Sprues.

Used to create a gateway for injection.

7014-JP



Matrix Pusher.

Used to facilitate putty application.

7010-JP



TCS Putty.

Used to create a matrix when injecting with the JP90 Handheld Injector. Can also be used to invest for clean injections.

Base & Catalyst
400gms ea.

4296-01

JP90 Starter Kits.

JP90 Starter Kits.

JP Starter Kit 7001-JPU

JP Starter Kit with Digital Furnace 7020-JPU

Includes TCS Digital Furnace



1 JP90
handheld
injector



JP bullets (small)
1 pk -
Unbreakable™
1 pk - iFlex™



JP bullets (large)
1 pk -
Unbreakable™
1 pk - iFlex™



JP flask



TCS putty
400 gms



2 waxed cases
for practice



1 matrix pusher
5 metal sprues


JP90 Bullets.

Refer to the table below for
bullet shades, sizes, and item
codes.

Unbreakable™ page 9

iFlex™ page 13

Karadent™ page 17

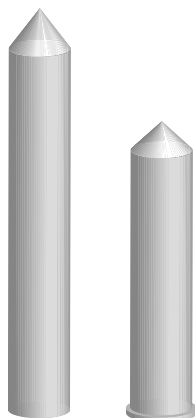

Digital Furnace & Accessories.



- » Simple and practical.
- » Melts thermoplastic resin.
- » Compact design.
- » Reliable / low maintenance.

110V 3512-01

220V 3512-02

SPECS.

H: 6" L: 1.5" W: 11"

1.3 lbs. (0.6 kilos).

2 year warranty

Made in USA



Melting Cylinder with Handle.

Fabricated from high precision stainless steel to withstand warpage under extreme heat.

3703-01



Knock Out Stand and Rod.

Steel tools used for purging the melting cylinder from the remaining cartridge and bronze plug after injection.

Knock Out Stand. 3705-01

Knock Out Rod. 3707-01



Bronze Plugs.

Thick bronze plug is used when melting medium and small cartridges in the furnace.

Thin bronze plug is used when melting large cartridges in the furnace.

Bronze Plug Thin. 3709-01

Bronze Plug Thick. 3712-01

Rental Equipment Program.

On the fence?

Rent a TCS injection unit for 20 days prior to investing.



Try Before You Buy.

Contact us to get started:

sales@tcsdentalinc.com

Waxes / Flasks.

Palate Wax.

Save time waxing with pre-formed palates created with the recommended thickness of a TCS palate.

5015-01 80/box



Lingual Wax.

Pre-formed lingual bars are ideal for clasps, flange, and lingual connector areas.

5020-01 1 lb



Sprue Wax.

Easy to work with, does not get brittle. Ideal diameter for spruing. 5 gauge.

5025-01 1 lb




Sprue Wax 9.5.

Connect your sprue wax to the 9.5 Sprue Wax only at the tip of the flask opening to avoid unwanted plaster sticking to it.

5027-01 5 pcs

5026-01 1 lb (~30 pcs)

Each piece yields ~10 injections.

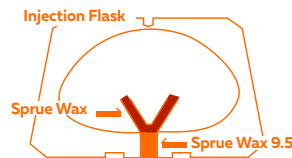

Injection Flask.

Milled from one piece of aluminum block. Extra smooth finish facilitates the de-flasking process. Includes a special knock-out button to simplify case removal and prolong the life of the flask.

3701-01



Pour Flask.

Used for acrylic pour technique and/or to duplicate models. Manufactured in light-weight aluminum. Easy to assemble, disassemble and clean.

3650-01



Nightguard Pour Flask.

Used for the fabrication of nightguards and/or duplicate models. Manufactured in light-weight aluminum. Easy to assemble, disassemble and clean.

3651-01



Brass Spruing Tools.

Used to sprue hydrocolloid paste for acrylic pour technique.

Medium (12mm dia) 3653M-01

Small (8.75mm dia) 3653S-01



Spatulas / Knives.

Plaster Spatula.

German steel blade with wooden handle. 4-1/4" flexible blade.

4230-01



Plaster Knife.

German steel blade with wooden handle. 2-1/4" blade.

4231-01



Fahnstock Knife and Spatula.

German steel blades with wooden handle for waxing.

4232-01



Wax Spatula.

German steel double ended spatula for waxing.

4233-01



Wax Spatula.

German steel double ended spatula for waxing.

4234-01



Wax Carver.

German steel double ended spatula and carver for waxing.

4235-01



TCS Carving Knife.

Excellent for removing extended fibers that sometimes appear when processing TCS partials. Also used to define the necklines of the TCS restorations when necessary.

4201-01



Replacement Blades.

Replacement blades for TCS Carving Knife.

4202-01 12 pcs



TCS Embroidered Apron.

Adjustable apron with pockets and embroidered TCS logo.

Apron01 1 pcs



Processing Tools.

TCS Cement.

To secure teeth in the stone prior to injection. Mild enough to keep teeth from bonding to the stone during injection.

3801-01 1/2 oz



Fusing Liquid.

For repairs and additions to existing TCS restorations. The Fusing Liquid is ONLY used when flasking and injecting are necessary.

3812-01 1/2 oz



Mold Release Spray.

Colorless and non-staining multi-purpose lubricating spray, used to ease the purging of the furnace accessories and TCS cartridges during injection process.

3814-01 12 oz



Hot Air Welder.

Hand torch in a slim design for minor repairs or additions to TCS restorations. Refillable with butane. (Butane not included).

3714-01



Twist Drills.

To create the diatoric retention in acrylic teeth necessary for fabricating TCS restorations.

005 5005-01 6 pcs

006 5006-01 6 pcs

009 5009-01 6 pcs

012 5012-01 6 pcs



EZ Flow.

Separating medium for processing thermoplastics. Eliminates stone sticking to the partial, leaving it shiny and clean.

4217-01 16 oz

4218-01 1 gallon



** Not recommended when hot flask injection is required, e.g. Karadent™.*

Empty Cartridges.

Empty cartridges are sold in packs of 50.



	25mm	28mm	22mm
Small	3003-01S		
Medium	3003-01M	3003-01MF	3003-01M22
Large	3003-01L	3003-01LF	3003-01L22

Finishing Tools.

Econo Cutter XL.

Use to cut off sprues from TCS restorations after injection.

4203-01



Econo Cutter Mandrel.

Mandrel required for the econo cutter.

4204-01



White Knock Down Wheel.

One of the first steps in the finishing stage, cuts excess bulk and trims partial to the desired length.

4211-01



White Knock Down Chuck.

Mandrel used with the White Knock Down Wheel and Rubberized Abrasives.

4212-01



Rubberized Abrasives.

Long-lasting wheel for smoothing the partial during the finishing stage.

Thin 4304-01

Thick 4306-01



Pink Grinding Wheels.

Excellent and economical product for the finishing stage.

4205-01 100/box



Brown Rubber Wheels.

Use prior to pumicing in order to obtain a high shine. It is crucial that the entire partial is rubberized and smoothed.

4503-01 100/box



Mini Fuzz Buff.

Use after the Brown Rubber Wheels to further smoothen partial prior to pumicing.

4222-01 12 pcs



Finishing Tools / Polishing Tools.

Brown Rubber Points.

Ideal to eliminate scratches and smoothen the surface.

Small 4513-01 12 pcs

Large 4514-01 12 pcs



Green Mounted Stones.

Used during the finishing stage. Cuts thermoplastics moderately, and leaves a smooth surface. Ideal for minor chairside adjustments.

4501-01 100/box



Carbide Bur.

Specifically designed to cut thermoplastics. Long-lasting with high cutting efficiency.

50251G-01



Pink Zirconia Fine Bur.

Use during the finishing stage to leave the surface smooth and prepped for the Brown Rubber Wheel. Also, ideal for chairside adjustments. Max RPM= 12,000. Use with a hand piece.

5035Z



Pink Zirconia Coarse Bur.

Use during the finishing stage to cut excess material and trim appliance to the desired length. Gross reduction. Max RPM= 16,000. Use with a hand piece.

5045Z



Polycril.

An advanced formula for pumicing acrylics and thermoplastics at high-speed.


4298-01 1 lb



Turbo Buff.

Pumicing buff with siliconized cotton cloth. Use for pumicing only.

4208-01

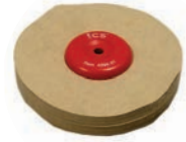

Rag Buff.

Stitched only in the center for the extra softness needed for polishing TCS restorations.

4207-01



Polishing Tools.

Mini Rag Buff.

Recommended to dentists for polishing touch-ups after minor adjustments on TCS appliances. Use with Green Hi-Shine polishing paste.

4221-01 12 pcs



Brown Tripoli.

Use after pumicing to prepare appliance for polishing. Effective for all TCS resins. Use with a dry Rag Buff.

4206-01 2 lbs



iShine.

Use after Brown Tripoli to polish surface of iFlex™ partials (and other polyolefin thermoplastics). Use with a dry Rag Buff.

4220-01 10.75 oz



Green Hi-Shine.

Use after Brown Tripoli to polish surface of Karadent™ and Unbreakable™ partials. Use with a dry Rag Buff.

4209-02 1 lb



Sample Model Unbreakable™.

An ideal marketing tool for your laboratory.

Upper Partial 4011-01U

Lower Partial 4011-01L



Sample Model: Karadent™.

An ideal marketing tool for your laboratory.

Upper Full Denture 4016-01U

Lower Partial Denture 4017-01L

Upper Clear Framework 4018-01U



Sample Model: iFlex™.

An ideal marketing tool for your laboratory.

Upper Partial 4009-01U

Lower Partial 4009-01L



Dr. Adjustment Kit.*

Ideal for quick and easy chairside adjustments.

4014-02



*Includes: 1 Green Point, 1 large Rubber Point, 1 small Rubber Point, 1 Pink Zirconia Fine Bur, 1 Carving Knife, 1 sample of Green Hi-Shine, 2 Mini Rag Buffs, 5 Mini Fuzz Buffs.



A clean appliance is vital to patients' health and the appearance and longevity of dental appliances. TCS Concentrated Dental Appliance Cleanser is a safe solution for quick and effective results.



TCS includes patient care instructions and a sample pack of denture cleaner with each cartridge.*

Include them with the finished case so that the patient can maintain a clean and vibrant appliance.

*USA only



Ideal for Dental Offices



TCS Concentrated Dental Appliance Cleaner.

Professional strength cleanser. Deep cleans, removes stains and odors. Ideal for most types of removable dental appliances.

4505-02 24 pcs
6 month supply

TCS Dental Appliance Cleaner.

TCS Dental Appliance Cleaner sample packets. Box of 75 packets (1-week supply ea.).

4505FR-SMP 75 pcs

TCS Dental Appliance Brush.

Resilient soft bristle brush to thoroughly clean dental appliances without scratching the surface.

4506FR-01



FREE DISPLAY WITH ORDER OF 9+ BOXES



Get More with TCS

- Technical Support.
- Equipment Rental.
- Training.
- Promotional Material.
- TCS Gives Back.

We Value Your Success

Quality tools and support help improve your laboratory's productivity and quality when fabricating a TCS Dental Appliance.

Connect with us to stay up to date on our latest information, videos and techniques.


Subscribe to our newsletter:

www.tcsdentalinc.com

Follow Us:



@tcsdentalinc

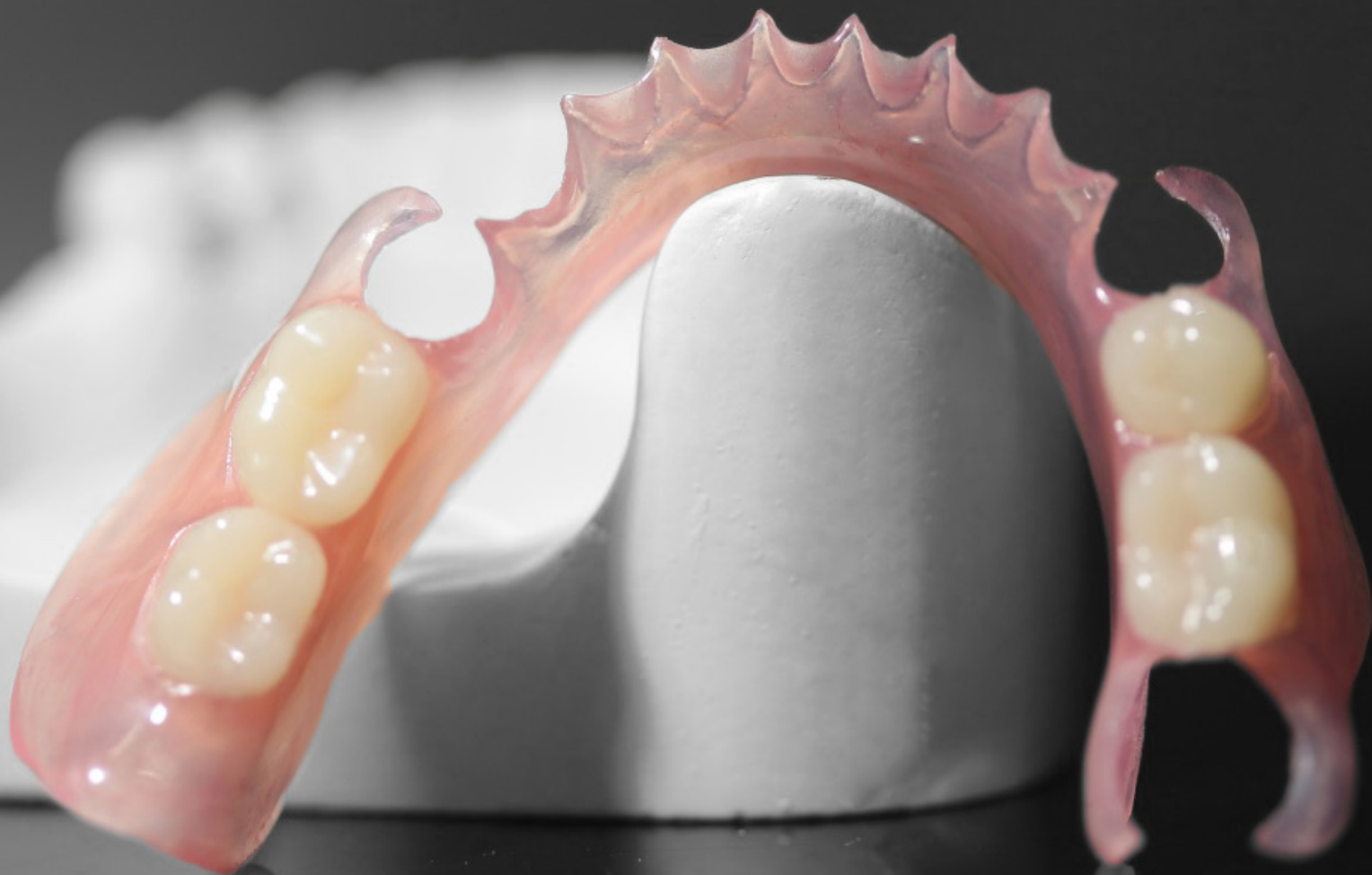

Technical Support.

How to find the most recent videos on our You Tube channel and register for TEXT Support.


YouTube Channel: *tcsdentalinc*

Subscribe to view short clips that will help with production and quality.

Playlists:

QUICK TIPS

1-3 min videos

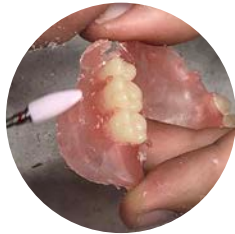

DESIGN LIBRARY

Design & blockout



SHORT VIDEOS

TCS Staff shares news



EQUIPMENT

Animated equipment demonstration



Text Support.

Register your phone for quick text support.

Text your name, city & state to (562) 241-1225.

Phone.

(866) 426-2970 | (562) 426-2970



TECHNICAL
QUESTIONS?

Text us!
(562) 241-1225

eLearning

Online training.

Learn how to fabricate TCS Flexible Partial dentures the TCS Way, wherever you are.

- Covers All the Basics**
Covers design, block out, setup, waxing, processing, finishing and polishing.
- Basic Training Kit**
Kits are available for USA customers only and includes all the basic tools needed to follow along with the course.
- Test your Knowledge**
Each chapter includes a quiz so you can test how much you've learned.
- Certificate of Completion**
Receive a digital certificate of completion upon completing the course.
- Earn CDT Credits**
The training course gets you 2 Scientific CDT credits!

Get Started Today
www.tcsdigitaltraining.com



Promote Your Removable Department.

Grow your removable department by utilizing these simple marketing tools.



**Flexible Partials
You Can Trust.**



Digital Media.

In today's digital age one of the most effective ways to promote your removable department is through online content. Enhance your website / social media platforms by embedding and sharing our digital marketing materials.

Contact us or visit our website for updated videos.

Ask about our banner / social media look book to utilize digital ads and more.

ASK

**about becoming
a TCS Certified
Laboratory**

JOIN • BE REFERRED • EARN

Contact us today for more information.

TCS Sales

sales@tcsdentalinc.com

Order Free Promotional Material



Patient Education Brochures.

Informative brochure for patients to learn about TCS appliances as a solution for missing teeth. Distribute to your accounts to display in their waiting area.

Size: 8.5" x 11".
Folded: 8.5" x 3.6".

A-8010-A1



Dr. Informational Booklet.

Includes everything a doctor needs to know about TCS appliances. E.g. impression taking, insertion, chairside adjustment techniques. A must have for dentists.

Size: 8" x 8".
(32 page booklet).

A021413



Dr. Flyer.

Convenient and informative pamphlet to distribute to your accounts. Includes information on all TCS resins & tools for their practice.

Size: 9" x 16".
Folded: 5.3" x 9".

122010



Your Name Goes Here.

Statement Stuffers.

Include these product stuffers in account statements. Includes space for you to add your logo and contact info.

Size: 8.5" x 3.5".

iFlex™ A020613
Karadent™ A040511
Unbreakable™ A033116

Feel free to contact us for promotional materials, digital material for websites / social media, and assistance with clinical studies.

CONTACT US TODAY.

Phone: 866-426-2970

Online: www.tcsdentalinc.com (Labs - Promotional Material).

Ships to USA only.

FREE OF CHARGE



Shade Guide.

Provide a shade guide to your accounts as a reference to match patient's natural gingiva to TCS resin shades.

Size: 3" x 2.5".

8015-02

together we smile

TCS GIVES BACK





Dominican Republic, 2019

TCS travels to remote areas across the globe to deliver partials and dentures to people in need. Our outreach program is possible thanks to our collaboration with participating local dentists and technicians. We are always looking for new opportunities.

Contact us for possible collaborations:

info@tcsdentalinc.com





2619 Lime Ave. Signal Hill, CA 90755, USA

(866) 426-2970 | (562) 426-2970

@**tcsdentalinc**

www.**tcsdentalinc**.com